
eppendorf

Move It[®] – Double Your Performance

Multi-channel Pipettes with Adjustable Tip Spacing

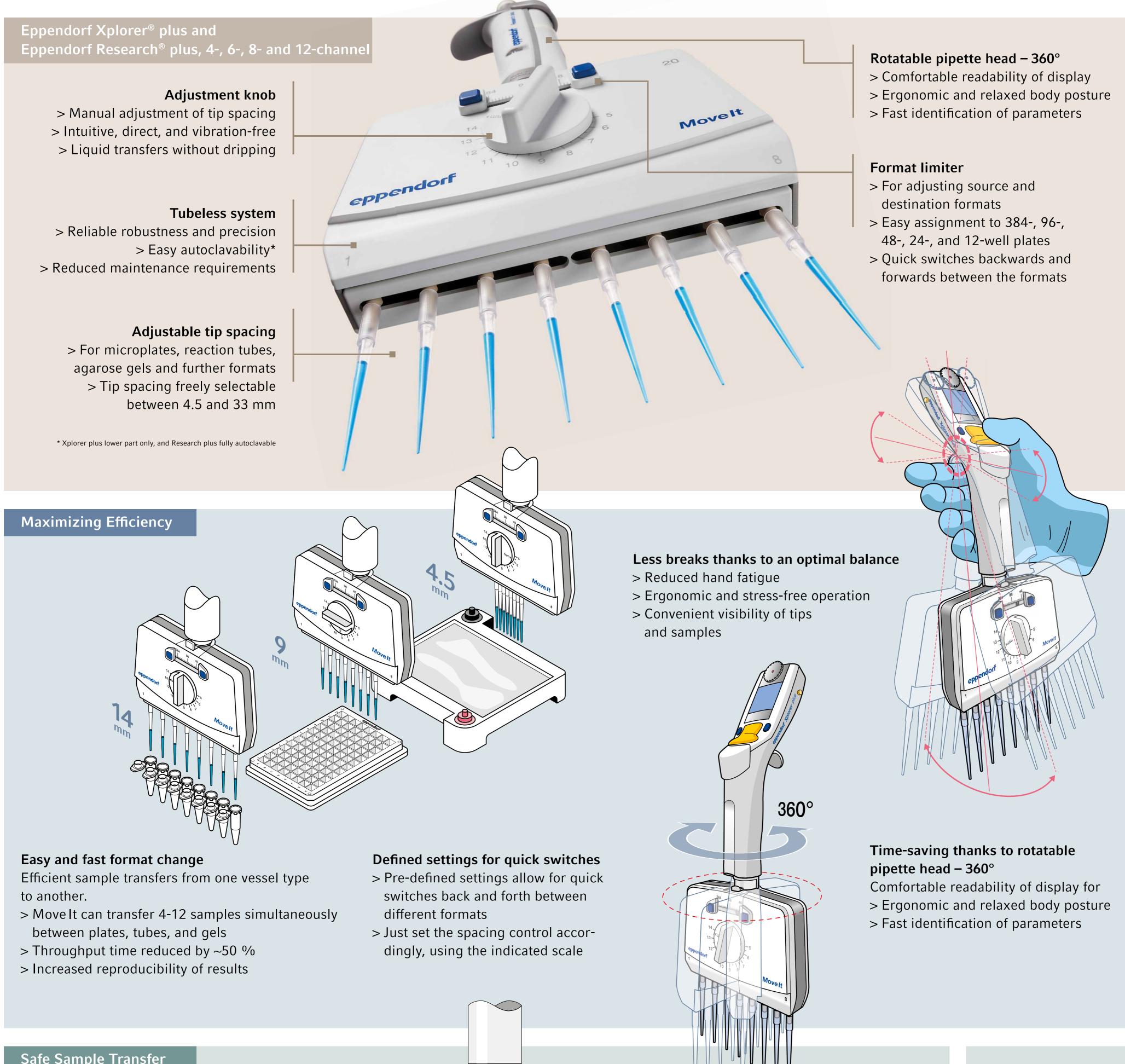

Safe Sample Transfer


Example applications

Tubeless system for reliability and reproducible results

Piston cylinder system directly connected to the nose, minimizing size of air-cushion.

- > Higher stability of air-cushion temperature and thus higher precision and reproducibility
- > Easy autoclavability* for reduced contamination risk
- > No cable tangle and tubing-porosity for consistently reliable robustness and precision
- > Less maintenance requirements

* Xplorer plus lower part only, and Research plus fully autoclavable

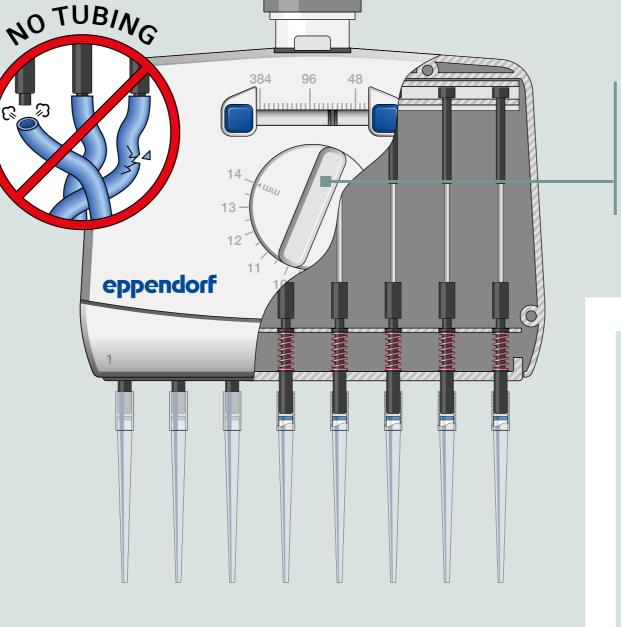

Adjustment knob

> Format changes without vibrations > Liquid transfer without dripping

> Blood analysis > FACS analysis > Cell Based Assays > Cell Culture > ELISA > PCR > Mouse genotyping

Move It selection guide

All models are available as electronic Xplorer plus and mechanical Research plus.

Type of tips	epT.I.P.S.®			epT.I.P.S.® 384	
No. of channels	4	6	8	8	12
Volume (µL)	300 / 1,200	300 / 1,200	300 / 1,200	20 / 100	20 / 100
Tip distance (mm)	9–33	9–20	9–14	4.5–14	4.5–9

www.eppendorf.com/move-it