thermoscientific

SPECIFICATION SHEET

Thermo Scientific Nunc EasYDish cell culture dishes

Novel ergonomic design improves ease of use

Thermo Scientific™ Nunc™ EasYDish™ cell culture dishes are designed to improve the handling, stacking, and transporting of cell cultures in your lab. The beveled grip ring (Figure 1) is easy to hold and stays secure in your gloved hand. The proprietary notch design (Figure 2) allows for optimal ventilation for cell growth. The raised lip (Figure 3) on the outer edge of the lid facilitates stable stacking of multiple dishes for secure transportation. The four orientation marks (Figure 4) on the bottom of the dish enable quickly locating cells under a microscope. Nunc EasYDish cell culture dishes come in four sizes (35 mm, 60 mm, 100 mm, 150 mm) along with a resealable package to maintain the cleanliness of opened packs.

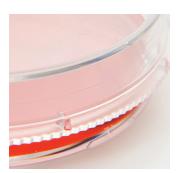


Figure 1. Beveled grip ring.

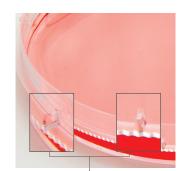


Figure 2. Proprietary notch design.



Figure 3. Raised outer edge for secure stacking.

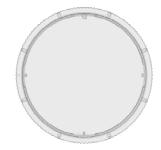


Figure 4. Orientation marks help you easily locate cells.



Feature	Benefits	
Beveled grip ring	Makes it easier to grasp and manage with a gloved hand	
Proprietary notch design	Maximizes ventilation while maintaining sterility	
	Minimizes friction between the lid and rim of the dish	
	Restricts the movement of the lid to help secure lid-to-dish fitting	
	Increases the physical strength of the dish	
Raised outer edge on dish lid	Allows multiple dishes to be stacked on top of each other with superior stability and security during transportation	
Orientation marks on dish bottom	Simplify locating cells under the microscope	
Writing area	Improves traceability of cell cultures	
Resealable packaging	Maintains cleanliness of opened packs	




thermoscientific

Product information	Specification	
Product testing	Free of RNase, DNase, and human DNA Noncytotoxic Nonpyrogenic Cell attachment lot testing	
Material	Polystyrene, USP Class VI	
Shelf life	5 years from manufacturing date	
Sterilization	Gamma irradiation to SAL 10 ⁻⁶	


Ordering information

Nunc EasYDish dishes with Thermo Scientific™ Nunclon™ Delta surface treatment			
Product	Cat. No.		
● 35 mm Dish, 8.8 cm², height 13 mm, sterile with lid, 10/pack, 500/case	150460		
● 60 mm Dish, 21.5 cm², height 16 mm, sterile with lid, 10/pack, 280/case	150462		
● 100 mm Dish, 56.7 cm², height 17 mm, sterile with lid, 10/pack, 150/case	150464		
● 100 mm Dish, 56.7 cm², height 21 mm, sterile with lid, 10/pack, 240/case	150466		
● 150 mm Dish, 145 cm², height 21 mm, sterile with lid, 10/pack, 80/case	150468		




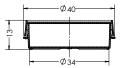


• 60 mm dish

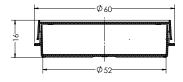


• 100 mm dish

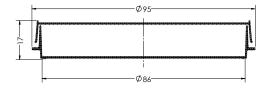

• 100 mm dish

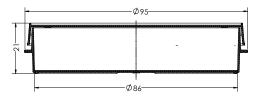

• 150 mm dish

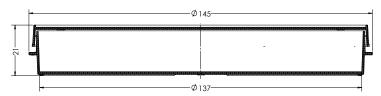

Explore EasYDish dishes at thermofisher.com/easydish

